

Quick Start Guide

Pairing Pure Charge&Go BCT IX to a PC



Bluetooth®
Connectivity
Transformed



Step 1

Turn your hearing aids off and back on again – this will put them in pairing mode for approximately 3 minutes.

Note: If your hearing aids are already paired to another device nearby, turn off the Bluetooth® on that other device during the pairing process.


Step 2

Open **Settings** by clicking the Windows icon in the task bar at the bottom of the screen, then clicking **Settings**.

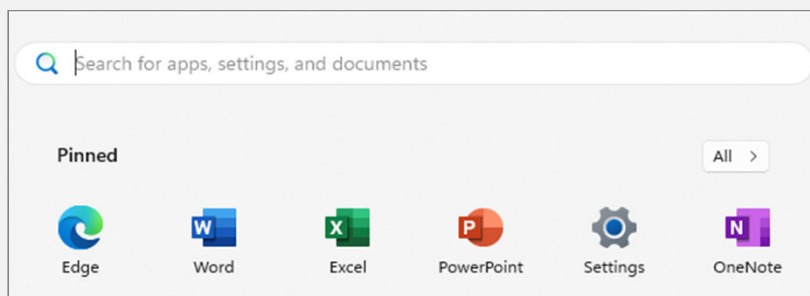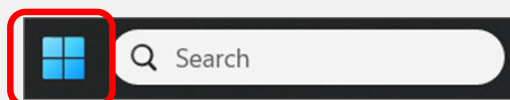


Step 3

On the left side of the Settings screen, select **Bluetooth & devices**.

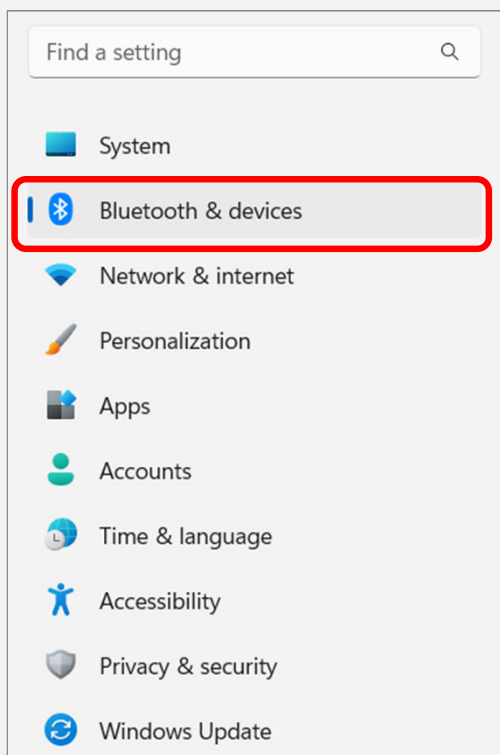

Step 4

Make sure the Bluetooth toggle is set to On and click **Add device**.

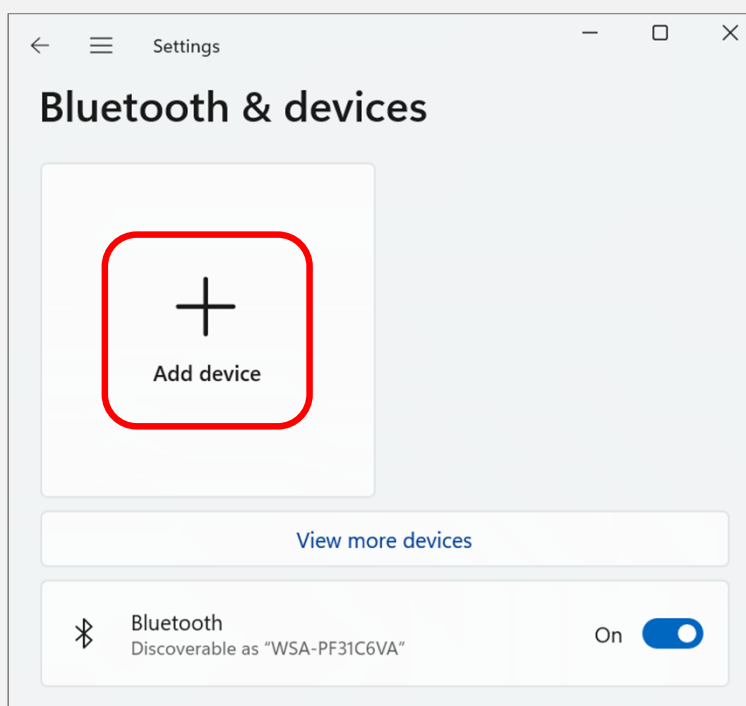

Step 5

Click **Bluetooth** device.

Add a device

Choose the kind of device that you want to add.



Bluetooth

Audio devices, mice, keyboards, phones, pens, controllers, and more



Wireless display or dock

Wireless monitors, TVs or PCs that use Miracast or wireless docks



Everything else

Xbox controllers with Wireless Adaptor, DLNA and more

Step 6

It may take a few moments for the hearing aids to be discovered. Click on the **hearing aid name** to complete the pairing process.

Add a device



Add a device

Make sure your device is turned on and discoverable. Select a device below to connect.



Anna's hearing aids

Audio

Step 7

Enjoy streaming from your computer to your Pure Charge&Go BCT IX hearing aids.

Add a device



Your device is ready to go!



Anna's hearing aids

Connected mic, audio

Disconnect



Bluetooth®
Connectivity
Transformed

Microsoft and Windows are trademarks of the Microsoft group of companies. The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by WSAUD A/S is under license.